
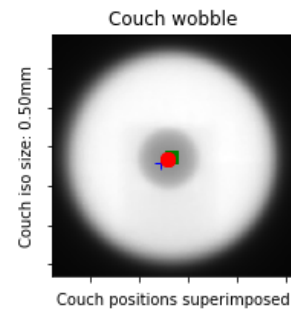
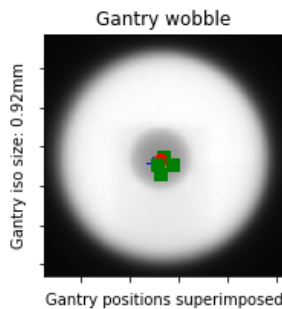
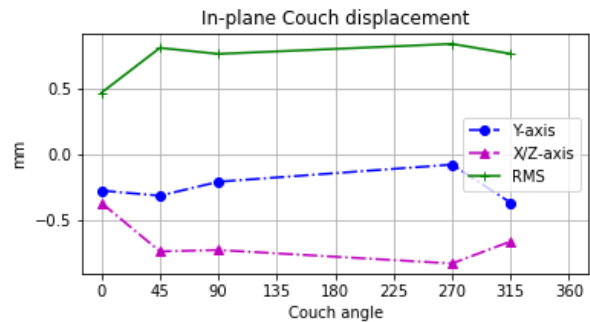
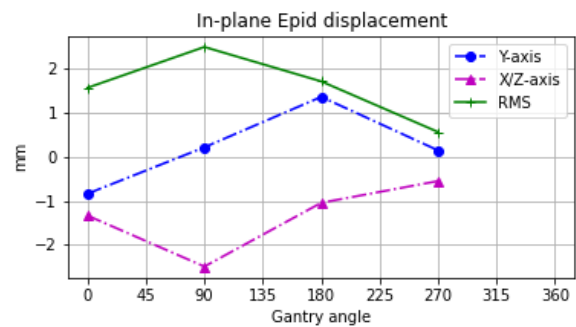



# Winston-Lutz Analysis

## Winston-Lutz Analysis

=====  
 Number of images: 8  
 Maximum 2D CAX->BB distance: 0.84mm  
 Median 2D CAX->BB distance: 0.72mm  
 Shift to iso: facing gantry, move BB: RIGHT 0.28mm; OUT 0.03mm; UP 0.20mm  
 Gantry 3D isocenter diameter: 0.92mm (4/8 images considered)  
 Maximum Gantry RMS deviation (mm): 0.67mm  
 Maximum EPID RMS deviation (mm): 2.50mm  
 Gantry+Collimator 3D isocenter diameter: 0.92mm (4/8 images considered)  
 Collimator 2D isocenter diameter: 0.00mm (1/8 images considered)  
 Maximum Collimator RMS deviation (mm): 0.00  
 Couch 2D isocenter diameter: 0.50mm (5/8 images considered)  
 Maximum Couch RMS deviation (mm): 0.84

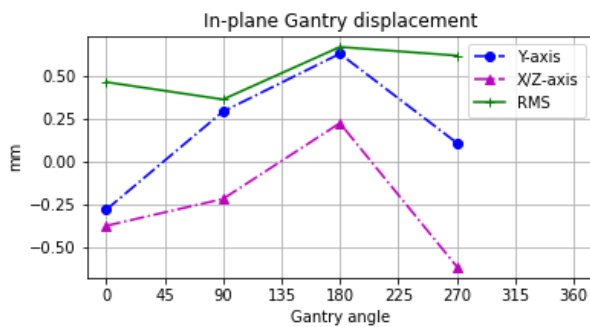

Gantry images

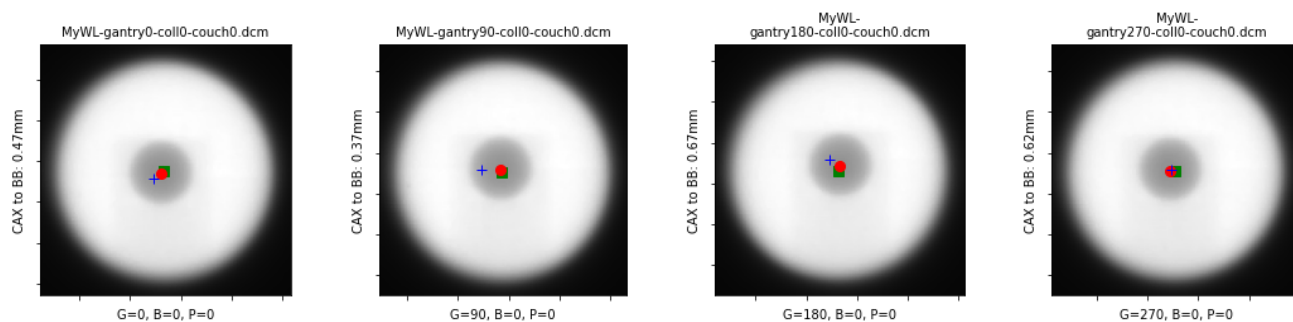

## Couch images

